
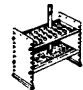














 <p>"E.R.L." HOLDERS for TEST TUBES, AMPOULES, BOTTLES, etc. In 8 sizes either in black synthetic rubber, or softer red organic rubber</p>	 <p>"ESCORUBBER" RUBBER and SILICONE RUBBER BUNGS In various qualities and many sizes</p>	 <p>"ESCORUBBER" ADJUSTABLE HOLDER for TEST TUBES & WATER BATH RACK Various sizes of holes available at 20/- each</p>
 <p>"ESCOPLASTIC" VINYL TUBING Non-pyrogenic and non-toxic for Blood Transfusion, sterilizable at 120 C., flexible, excellent for general laboratory and hospital work, milk, various chemicals, hot water, etc.</p>	 <p>"ESCORUBBER" SILICONE RUBBER TUBING Non-toxic, translucent, will withstand temperatures up to 200, 250°C.</p>	
 <p>"ESCORUBBER" RUBBER and SILICONE RUBBER BOTTLE CAP LINERS Made in many standard sizes or to specification</p>	 <p>"E.R.L." RUBBER FILTERING CONES for BUCHNER FLASKS In sets of 6 at 22s. 6d. per set</p>	 <p>"ESCORUBBER" RUBBER and SILICONE RUBBER SHEETS Made in various qualities and many sizes</p>
 <p>"E.R.L." SILICONE RUBBER COVERED MAGNETS for MAGNETIC STIRRING 1" at 7s. 6d. each 1½" at 9s. 6d. each 2" at 12s. 6d. each</p>	 <p>"ESCORUBBER" SR 300 SILICONE RUBBER MOULDING PASTE</p>	 <p>"E.R.L." RUBBER PIPETTE TEATS Extra heavy, in 7 colours for identification purposes approx. 2½ c.c. and 5 c.c. also approx. 7½ c.c. & 10 c.c. (red only)</p>
 <p>"ESCORUBBER" R 33 RUBBER TUBING High quality for various laboratory and hospital applications</p>	 <p>"ESCORUBBER" CLOSURES for BLOOD and SALINE SOLUTION BOTTLES Ref. No. P2—made in high quality rubber, and used extensively by many well-known hospitals</p>	 <p>"ESCOPLASTIC" PLASTIC PIPETTE TEATS With remarkable resistance to many chemicals In 7 colours for identification purposes at 45s. per gross</p>
 <p>"ESCOPLASTIC" POLYTHENE DISPOSABLE GLOVES 1000 at 16s. 8d. per 100 plus 1s. 4d. P.T.</p>	 <p>"ESCOPLASTIC" PLASTIC DISPOSABLE PETRI DISHES</p>	 <p>STAINLESS STEEL GRADUATED MEASURES and OTHER LABORATORY and HOSPITAL EQUIPMENT</p>
<p align="center">ESCO (RUBBER) LTD. Manufacturers of Rubber, Silicone Rubber and Plastic Products SEAL WORKS · SEAL STREET · LONDON · E.8</p>		

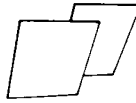6

easy steps


to the rapid, accurate diagnosis of pregnancy




Filter a small quantity of the urine to be tested.




Transfer the appropriate amount of buffer solution to vial of sensitized erythrocytes.

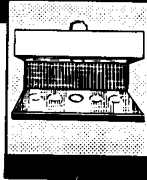

Break off neck of Pregnosticon ampoule and place in mirrored rack.

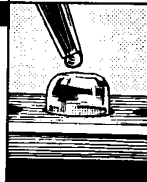


Transfer 0.1 ml. of the filtered urine to the ampoule.

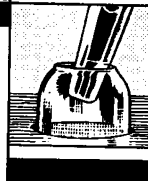


Add 0.4 ml. of the erythrocyte suspension. Shake and allow to settle free from vibration and heat.

READ RESULT 2 HOURS LATER

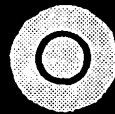

PREGNANT

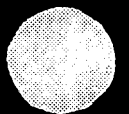

NON-PREGNANT

PREGNOSTICON

Immuno-diagnostic Pregnancy Test



Organon Laboratories Ltd., Crown House, London Road, Morden, Surrey

for proven accuracy, precision and reliability —



SGOT, SGPT

DADE TRANSAMINASE REAGENTS are used in a simple, direct, and proven procedure for colorimetrically measuring the levels of serum glutamic oxalacetic transaminase (SGOT) and serum glutamic pyruvic transaminase (SGPT). Each set contains stable, balanced reagents for 100 determinations.

DADE TRANSAMINASE REAGENTS are standardized to yield the same values as obtained by the U. V. reference method of Karmen and Wroblewski. You prepare only one calibration curve and verify it with Enza-trol[®], the known enzyme test control serum.

LDH REAGENTS

DADE LACTIC DEHYDROGENASE (LDH) REAGENTS are used in the determination of serum or body fluid LDH levels. Elevations of LDH are often demonstrated in cases of infarctions or malignancies. The test is conveniently performed in a vial containing DPNH.

LDH and SGOT tests, when run in parallel, provide valuable additional information as an aid in the diagnosis of myocardial infarctions.

HBD REAGENTS

DADE HYDROXYBUTYRATE DEHYDROGENASE (HBD) REAGENTS, for the quantitative colorimetric measurement of HBD activity, have a diagnostic advantage over determinations of total serum LDH and SGOT in cases of myocardial infarction. The peak activity of HBD is obtained at the same time as LDH but remains elevated for a longer period than either LDH or SGOT. HBD has been reported to be specific for indicating heart tissue damage.

HBD determinations, when performed in conjunction with the LDH test, can differentiate between liver damage and heart tissue damage.

ENZA-TROL[®]

DADE ENZA-TROL[®] fulfills the need for a known enzyme test control serum. It provides you with a reliable check on substrates, reagents, instrument variation and technic. Used exactly as serum, Enza-trol[®] contains serum protein and known values of enzymes normally found in human serum. Because most enzymes are unstable in solution, Enza-trol[®] is supplied in the dried form. After reconstitution, stability is the same as enzymes in human serum.

DADE

AUTHORISED DISTRIBUTORS IN THE UNITED KINGDOM

ARMOUR PHARMACEUTICAL COMPANY LTD / Hampden Park / EASTBOURNE / Sussex / Eastbourne 51111

ELECTRICALS (LABORATORY) LTD / 'Hedgefield House' / Blaydon-on-Tyne / Co DURHAM / Ryton 2381

ELECTRICALS (LABORATORY) LTD / 3A St Vincent Street / EDINBURGH 3 / CAL 2188

GEORGE T. GURR LTD / 136-140 New Kings Road / LONDON SW6 / REN 5482

ARNOLD R. HORWELL LTD / 17 Cricklewood Broadway / LONDON NW2 / GLA 8777 & 8631



KERNICK & SON LTD / Monthermer Road / CARDIFF / South Wales / Cardiff 26128

ORME SCIENTIFIC LTD / Bentley Avenue / Middleton / MANCHESTER / Middleton 2648

E. J. BROWNE, M.P.S.I. / 6 The Crescent / Sallynoggin Dunloaghair / DUBLIN / EIRE / Dublin 87020

— you can depend on Dade

Reagents for
coagulation tests

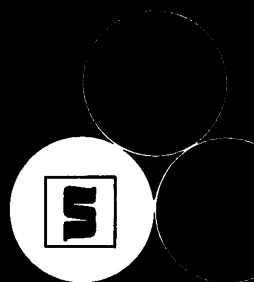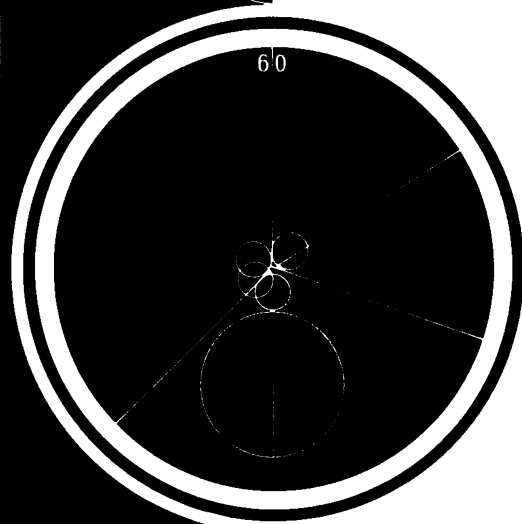

**Deprothrombinised
Ox Plasma**
(as used for OWBENS
P&P method)

Price: £6 per 100 vials
each containing
2 ml Ox Plasma
£4 per 50 vials each
containing 2 ml
Ox Plasma
£2 per 25 vials each
containing 2 ml
Ox Plasma

**Rabbit Brain
Thromboplastin
(Acetone Dried)**

Stability proven.
Results Reproducible.
Price: £1.15.0d. per box
of six vials—each vial
sufficient for 5ml
of Thromboplastin.

**Rabbit Brain
Thromboplastin
(Saline Extract—
Phenolized)**

Price: £1.15.0d per box
of 6x5 ml bottles.

**Bell & Alton
Platelet Substitute**
(For Thromboplastin
Generation Tests)

Price: 18/- per box
of 6x5 ml bottles.




STAYNE
STAYNE
STAYNE
STAYNE
STAYNE
STAYNE


Stayne Laboratories Ltd
Hillbottom Road, Sands,
High Wycombe, Bucks.
Tel. High Wycombe 71, 1

TUBERCULIN TINE TEST


(Rosenthal) Lederle

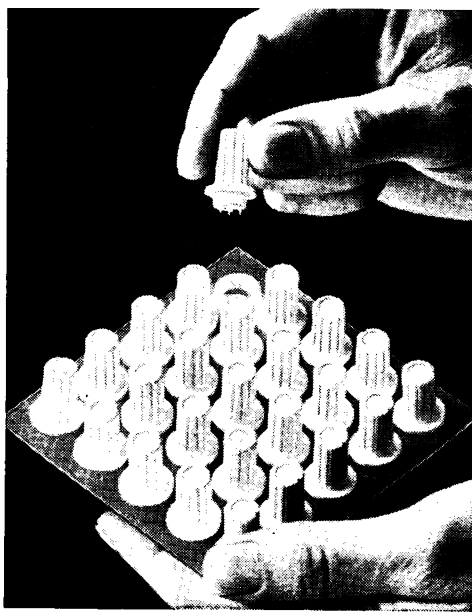

A PRACTICAL ADVANCE IN TUBERCULIN TESTING

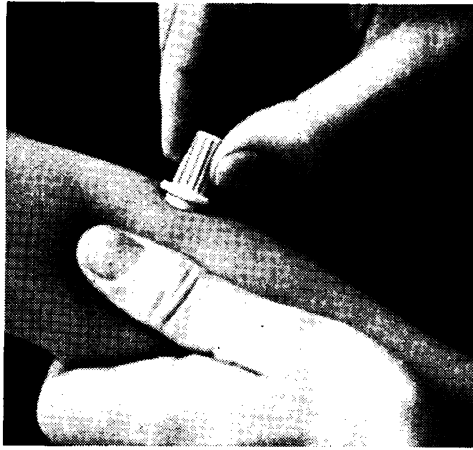

- **Simple** - requires no additional equipment
- **Accurate** - standardised against intradermal Mantoux test (5 T.U.)
- **Convenient** - small, portable and ready for immediate use.
- **Time-Saving** - each test takes only seconds.

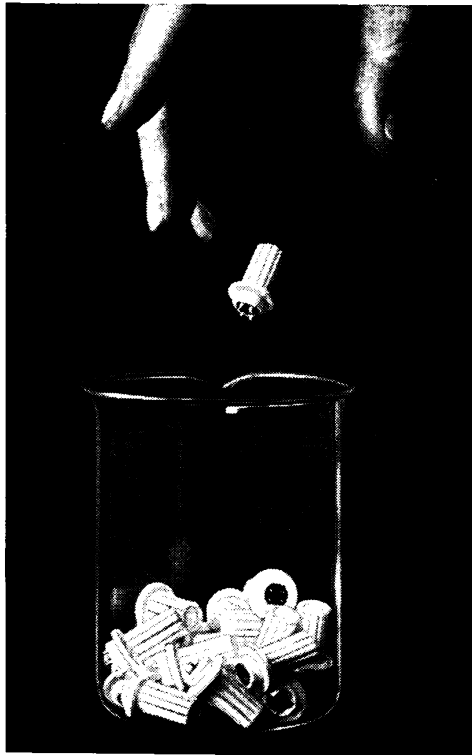

- **Disposable** - no soiled apparatus no clearing up.


box of 25 tests

LEDERLE LABORATORIES

a division of
Cyanamid of Great Britain Limited



CYANAMID

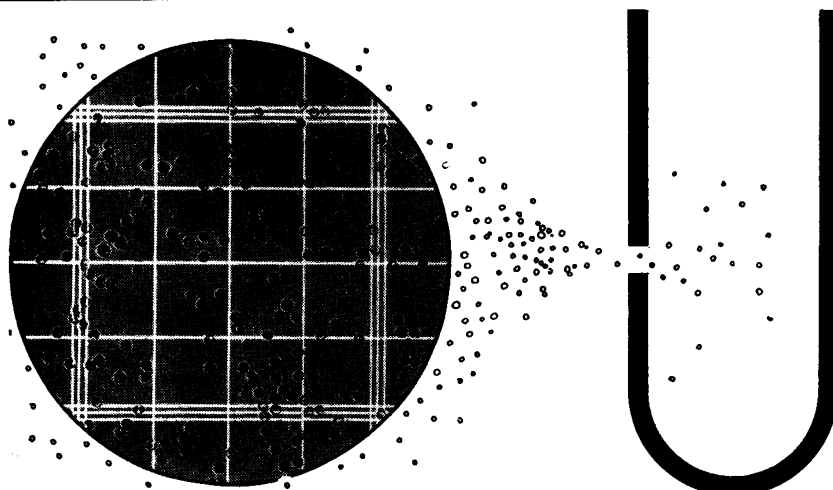

Red blood cells
Platelets
Algae
Tissue culture cells
Bacteria
Plant cells
White blood cells
Plankton
Protozoa
Bone marrow cells
Mitochondria
Haemagglutins

Whether your interests are solely the routine counting of red and white cells or entail size distribution studies of bacteria, platelets, tissue culture cells, spermatozoa or other biological material, the modern approach is the Coulter Counter. There is a model to suit your specific needs and priced to suit your finances . . . Five Thousand cells counted and sized every second with a simplicity and accuracy second to none.


**COULTER
COUNTS
FASTER**



TRADE MARK

COULTER ELECTRONICS LIMITED

Ashwell Street, St. Albans, Herts. St. Albans 50828

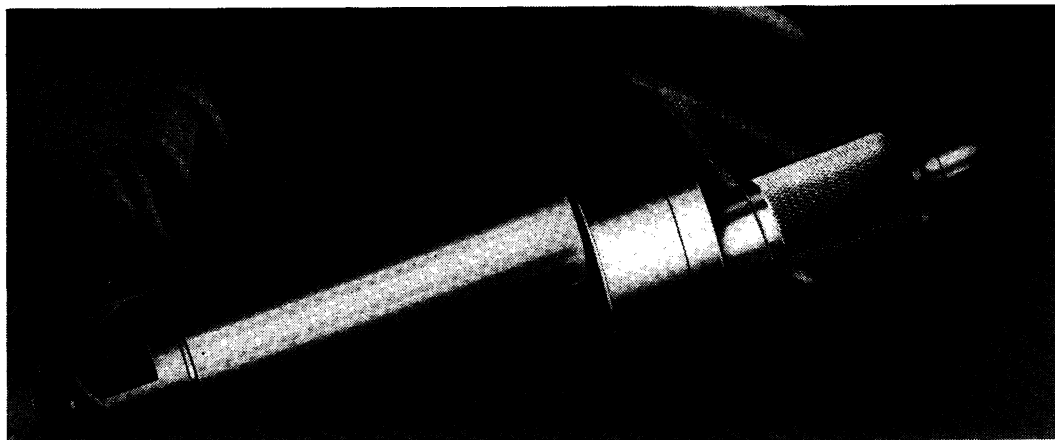

a saw that doesn't rotate!

The latest and most refined version of bone saw, the Desoutter Pneumatic oscillating bone saw (weight 1½ lbs.) has a high speed movement of just a few degrees each way, will cut only firm substances on positive contact. Greater safety is ensured and spraying of

liquid and fragments is eliminated. Speed is precisely controlled by light action trigger and there is ample power for cutting bone grafts. This instrument represents a really significant advance in orthopaedic surgery. Descriptive literature gladly sent on request.

Desoutter

PNEUMATIC OSCILLATING BONE SAW
COMPLETE WITH ACCESSORIES